

IN THE WORLD OF PHARMA PROCESS SOLUTIONS.

WE GO RIGHT AROUND THE WORLD!



We are anish pharma - a globally renowned name in Pharma Process Solutions. In three decades, we've reinvented ourselves and evolved with the rapidly-changing times. Today, we are an end-to-end process solutions provider for OSD, Semi-solid & Liquid. We cater to the Pharma, Nutra and Food industries with installations across continents.



35+
Years




30+
Countries



1,550+
Installations



1,00,000+ Sq Ft
Mfg. Area

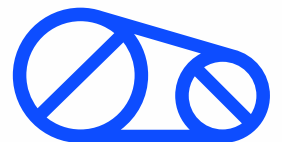

anish

anish granulation line



1 L to 2500 L scalable range



High yield and less product loss due to closed transfer by dust free charging and discharging



Through-the-wall & Stand-alone installation options



Independent operating interface with integrated PC controls for automated processes and short set-up times



Convertible combo processor: Drying, Top spray & Bottom spray options



IPC Bin Lifter, Tipper & Multi-Purpose Column options for product charging & discharging



Wet & Dry Co-Mills, Vibratory & Rotary Sifter options



Octagonal, Double-Cone & Conta-Piller Blender options



anish mixer granulator



1 L to 2500 L complete scalable range



Inline under-driven & over-driven
co-mill at discharge



Dust free charging by vacuum
& air purged degassing filter



Direct & belt driven
impeller technology



Advanced mechanical seals
with air purging arrangement



Binder addition with spray nozzle
& peristaltic pump



Jacketed bowl option for
hot melt granulation



Through-the-wall &
Stand-alone design options



anish extruder & spheroniser

Extruder



Output range from
0.5 kg/hr to 400 kg/hr



Single & twin screw extruder options



Easily convertible in
Dome / Radial / Axial extruder



Jacketed design for heating
& cooling provision



Spheroniser



Available in disc diameter
200 mm to 900 mm



Chequered plate options
from 1 mm to 8 mm pitch



Plate lifting mechanism with
bucket for easy dismantling



Jacketed design for heating & cooling



GAMP 5
Category 4 Compliant



anish fluid bed



AFB - C



ATSG - C




APCG - C




Drying



Top Spray

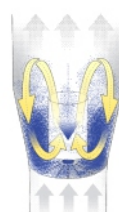

Bottom Spray

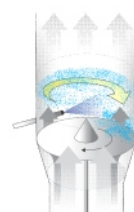

Tangential Spray



1 L to 2500 L complete scalable range



Cylindrical & conical expansion chamber design options



Convertible combo processor:
Drying, Top spray, Bottom spray
and Roto granulator options



Single shaking & twin shaking
filter bag and air purge socks
& bonnet filter options



Top spray granulation-1, 3, 6 & 9 head spray
nozzles. Bottom spray coating - Standard,
high speed & light design spray nozzles.



Precise dew point control
(up to -40°C) to process
moisture sensitive products



Through-the-wall &
Stand-alone design options



Product charging and discharging
by vacuum transfer and tipper options



anish tablet coater



1L to 2000L scalable range



Ideal for film coating & sugar coating on tablets



10% to 100% batch occupancy



Dust free charging & ensures 100% discharge



Contactless product temperature measurement



High efficiency manual & robotic coating spray arm options



GAMP 5
Category 4 Compliant



anish containment system



Fully closed ISO Class 8 environment with
CLASS 2 Leak Tightness classification-
Tested according to ISO 10648-2



IPC SCADA or HMI-based operation
for user-friendly control



Suitable for Containment up to OEB 5



Comprehensive High Containment
Accessories, including Transfer
Valves and Bags



Tailored for versatile applications:
sampling, dispensing, granulation,
tablet coating, etc.



Rapid Transfer Technology for
efficient Inter-Isolator Connectivity



Hassle-free cleaning with
an integrated Wash in Place System



Options for VHP or Nitrogen
Inertization for added versatility



anish liquid manufacturing plant



10L to 10000L scalable range



Inline, inbuilt top & inbuilt bottom type emulsifier options



Total process plant including premix, active phase mixing and liquid manufacturing technology



Load cell for precise weighing



Top entry & Bottom entry stirrer options



Completely automated process with vacuum product transfer



GAMP 5
Category 4 Compliant



anish ointment manufacturing plant



10L to 5000L scalable range



Inline or inbuilt type of high shear emulsifier options



Total process plant including premix, active phase mixing and ointment manufacturing technology



Load cell for precise weighing



Top entry & Bottom entry anchor options



Completely automated process with vacuum product transfer



GAMP 5
Category 4 Compliant



anish mixing vessel



Capacity range 5L to 500L for mobile vessel and 10L to 3000L integrated vessels.



Bottom-mounted magnetic mixers provide efficient mixing with low, medium and high shear options.



Flush bottom outlet valve ensures optimized drainage with customizable sampling and sterilization options



Sterile flange connections on the top dish and side ports enhance CIP/SIP compatibility.



Integrated sanitary ports support pH, temperature, conductivity and DO monitoring & control.



Load cells and level sensors ensure precise weight and level measurement.



Interconnected product piping through sampling systems and pumps with a filtration system ensures efficient filtration for integrated vessels.



global presence



NORTH AMERICA - USA • Mexico

SOUTH AMERICA - Brazil

AFRICA - Algeria • Nigeria • Ghana • Tanzania • Kenya • Tunisia • Uganda • Sudan

EUROPE - Russia • Cyprus

ASIA - India • China • Philippines • Vietnam • Malaysia • Iran • Pakistan • Bangladesh, UAE • Indonesia • Saudi Arabia • Jordan • Syria • Iraq • Thailand • Sri Lanka

AUSTRALIA - Australia • New Zealand



anish pharma equip pvt ltd

HEAD OFFICE:

G-68, M.I.D.C, Ambad, Nashik,
Maharashtra, India - 422010.

W www.anishpharma.com

SALES OFFICE:

B-613, Lodha Supremus II, Wagle Estate,
Thane (W), Mumbai, Maharashtra, India - 400604.

E marketing@anishpharma.com